HD Balance





The tilt-in-space wheelchair HD Balance is avaliable in two different wheel sizes 16" pushing wheelchair 24" self-propelled wheelchair

four seat widths 38, 42, 46 and 50 cm. The wheelchair offers several accessories and can be costumized to the wishes of the client.

Features:

Handlebar:

height adjustable

Back:

Flexi back with velcro straps

Backrest:

angle adjustable

Seat measurements:

Seat width & seat length adjustable

Armret:

height adjustable, detachable, upholstered

Legrest:

angle adjustable

Footrest:

divided, angle adjustable

Brakes:


Attendant brake


Other:

calf support



Flexi Back

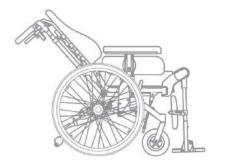

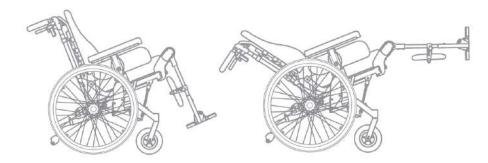

Leg support with extremely good agle adjustment

The leg supports can be angled between 90 and 180 degrees by default, without the casters need to be replaced with smaller ones. It is also very easy to remove or swing away the leg support unit. As an option there are leg supports available in different fixed angles.

Robust and flexible armrest

The armrest is mounted in one anchor point which makes it easy to take off and put in place. Setting options in depth and height is generous. The armrest also features a mounting rail where various accessory brackets can be mounted, for example table mounts.

A wide and flexible range of accessorie

A variety of accessories are available such as alternative leg and calf support, trunk support, shoulder support, tables etc.

The wheelchair has a clean underbody and high ground clearance



Good Handling

Extremely short distance between the caster & wheel gives very good handling for both users & carers. The driving wheels are also fully in line with the castors allowing the chair to run easily on narrow ramps.

Seat tilting pivot is located near user CG

Provides better control of tilting as it requires less force to tilt, especially valuable when tilting heavy users. The tilting pivot position offers the opportunity for expanded tilting in special cases.

High patient security

Scales for different settings

The good seating ergonomics is combined with a number of scales to make it easier for carers and relatives to ensure that the user is sitting in the right position. The communication between prescribers and carer / user / family is greatly simplified by this.

Control Placement

Tilt and back angle controls are positioned so that the carer can maintain his/her grip around the push bar when using these functions. It gives greater control of the movement!

Flexi Back with velcro-straps

Flexi back with adjustable Velcro-straps is standard for HD Balance. As an option it is possible to choose increased length and/or width.

Integratred transport mounts

Telephone: +49 341 2612498 Online: www.cs-humanmobility.com

| Measures in [mm] if not stated | balance 38 16" 24" | balance 42 16" 24" | balance 46 16" 24" | balance 50 16" 24" |
|--|-----------------------|-----------------------|-----------------------|-----------------------|
| | | | | |
| Seat width | 380 | 420 | 460 | 500 |
| Seat depth standard frame, with cushion | 430-520 | 430-520 | 430-520 | 430-520 |
| Seat depth standard frame, without cushion | 380-470 | 380-470 | 380-470 | 380-470 |
| Seat depth extended frame, with cushion | 480-570 | 480-570 | 480-570 | 480-570 |
| Seat depth extended frame, without cushion | 430-520 | 430-520 | 430-520 | 430-520 |
| Distance back-knee joint, standard frame | 560 | 560 | 560 | 560 |
| Distance back-knee joint, extended frame | 610 | 610 | 610 | 610 |
| Seat height, without cushion | 450 | 450 | 450 | 450 |
| Seat heght with original cushion, rear end | 510 | 510 | 510 | 510 |
| Seat heght with original cushion, front end | 540 | 540 | 540 | 540 |
| Seat height, without cushion, increased height | 500 | 500 | 500 | 500 |
| Seat height, without cushion, lowered height | 420 | 420 | 420 | 420 |
| Seat height, without cushion, 20" / 5" wheels | n/a 410 | n/a 410 | n/a 410 | n/a 410 |
| Back height standard, seat plate-top of back frame | 550-600 | 550-600 | 550-600 | 550-600 |
| Back height extended, seat plate-top of back frame | 650-700 | 650-700 | 650-700 | 650-700 |
| Back width | 380 | 420 | 460 | 500 |
| Armrest height standard, from seat plate | 255-345 | 255-345 | 255-345 | 255-345 |
| Armrest height low, from seat plate | 195-255 | 195-255 | 195-255 | 195-255 |
| Armrest depth, from back frame reference point | 340-430 | 340-430 | 340-430 | 340-430 |
| Armrest length | 370 | 370 | 370 | 370 |
| Armrest width | 70 | 70 | 70 | 70 |
| Leg support angle [º] | 90-180 | 90-180 | 90-180 | 90-180 |
| Footrest height (seat plate - foot plate) | 230-430 | 230-430 | 230-430 | 230-430 |
| Footrest angle (legrest tube-footrest) [°] | 90-110 | 90-110 | 90-110 | 90-110 |
| Footrest, depth | 210 | 210 | 210 | 210 |
| Footrest rear end - knee joint (centre) | 70-110 | 70-110 | 70-110 | 70-110 |
| Calf support height | 210 | 210 | 210 | 210 |
| Calf support depth, to centre of knee joint | 60-90 | 60-90 | 60-90 | 60-90 |
| Push bar height, standard, 3-positions | 1050-1200 | 1050-1200 | 1050-1200 | 1050-1200 |
| Total width | 630 640 | 670 680 | 710 720 | 750 760 |
| Total length, with legrests | 980 | 980 | 980 | 980 |
| Total length, without legrests | 760 | 760 | 760 | 760 |
| Total height, back in 90° position | 1050 | 1050 | 1050 | 1050 |
| Total height, no cushion & forward tilted back frame | 750 | 750 | 750 | 750 |
| Total weight [kg] incl. original cushions | 29,8 29,9 | 30,1 30,2 | 30,4 30,5 | 30,7 30,8 |
| Turning space | 1340 | 1340 | 1340 | 1340 |
| Seat tilt angle [°], (-5° possible) | 0 - 20 | 0 - 20 | 0 - 20 | 0 - 20 |
| Seat tilt angle, extended seat tilt [º] (-5º possible) | 0 - 30 | 0 - 30 | 0 - 30 | 0 - 30 |
| Back angle [º] | 90-120 | 90-120 | 90-120 | 90-120 |
| TANSPORTATION MEASURES (Collapsed, without; w | | | | |
| Width | 630 | 670 | 710 | 750 |
| Length | 750 | 750 | 750 | 750 |
| Height | 720 | 720 | 720 | 720 |
| Weight [kg] | 16,5 | 16,7 | 16,9 | 17,1 |

Telefon: +49 341 2612498 Online: www.cs-humanmobility.com